

Delta Plus Mattress

A foam support surface aids in the prevention and management of pressure injury. Pressure, shear, friction and moisture are all contributing factors effecting skin integrity.




Specifications


Product Code	DP19888Z1593
Length	1980mm
Width	880mm
Depth	180mm
Maximum Weight	250kg
Mattress Weight	17kg
Cover Warranty	3 years
Foam Warranty	7 years

Foam Specifications


- Delta Plus Mattress has been designed using three layers of hospital grade foam. Density, hardness, and depth play a major role in its construction.
- Top surface under the torso has been castellated increasing surface contact area and potentially reducing risk factor while the heel section has a gel memory foam insert.
- 3 layered foam
 - 60mm Top layer 39/140 Castellated
 - 60mm Top layer 55/50 Heel
 - 70mm 2nd layer 39/170
 - 50mm Base layer 36/200

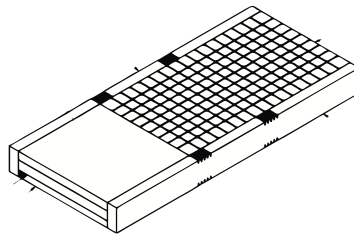

Options available:

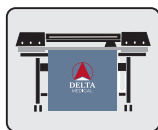

Custom printing

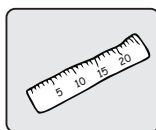

Custom sizing



Tracking



Product Code:
DP19888ZE1593

- Mental Health Mattress option also available: • **Product Code: DP19888MH1593**
- Fully welded seams
 - No zips
 - No handles
 - With air vents



All seams are welded ensuring it is impervious to fluid penetration.

Cover Specifications

Top cover

- Vapour Permeable
- Multi directional elasticity
- Waterproof
- Antibacterial
- BS 6807 ignition source 5
- BS 7175 ignition source 5
- Fully welded with zip (4 sides) and waterfall flaps
- 3 year guarantee on fabric if maintained correctly
- Care instructions on base cover

Base Cover

- Grey vinyl non slip for durability and longevity
- Carry handles

Standard Colour



Carrflex Extreme Royal Blue



DELTA
MEDICAL