

Previcare Mattress

Delta Medical's Previcare Mattress consists of 3 layers of foam offering support from the base layer to superior comfort from the castellated and gel memory top layer foam.




Specifications


Product Code	PC19886Z1660
Length	1980mm
Width	860mm
Depth	150mm
Maximum Weight	200kg
Mattress Weight	14kg
Cover Warranty	3 years
Foam Warranty	7 years

Foam Specifications


- Delta Medical Previcare Mattress has been designed using three layers of hospital grade foam. Density, hardness, and depth play a major role in its construction.
- Top surface under the torso has been corrugated increasing surface contact area and potentially reducing risk factor while the heel section has a gel memory foam insert.
- 3 layered foam
 - 50mm Top layer 39/115 Corrugated
 - 50mm Top layer 55/50 Heel
 - 70mm 2nd layer 39/140
 - 30mm Base layer 36/200

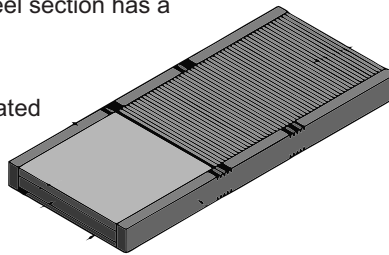

Options available:

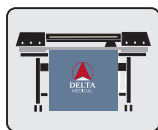

Custom printing

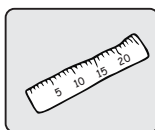

Custom sizing



Tracking



Product Code:
PC19886ZE1660

- Mental Health Mattress option also available:
- **Product Code: PC19886MH1660**
 - Fully welded seams
 - No zips
 - No handles
 - With air vents

Cover Specifications

Top cover

- Vapour Permeable
- Multi directional elasticity
- Waterproof
- Antibacterial
- BS 6807 ignition source 5
- BS 7175 ignition source 5
- Fully welded with zip (4 sides) and waterfall flaps
- 3 year guarantee on fabric if maintained correctly
- Care instructions on base cover

Base Cover

- Grey vinyl non slip for durability and longevity
- Carry handles

Standard Colour



Carrflex Extreme Sapphire Blue



All seams are welded ensuring it is impervious to fluid penetration.



DELTA
MEDICAL