

# Isolation Gown

A disposable non-woven, lint free, polyethylene coated fabric. This fabric offers maximum protection and is recommended for extended wear situations involving high fluid contact.

## KEY SPECIFICATION

Product Code 7056WCBM

Sizing Medium

Quantity 50pcs / carton

Material PP/PE 40GSM

ARTG# 153273

## STRUCTURE

This gown consists of two layers of non-woven fabric. Outside layer is polyethelene while the inner layer is spunbonded polyproplyene.

### SPUNBOND

Spunbond is a single layer non-woven fabric which is broadly defined as a web structure consisting of long continuous strand of filaments offering durability and strength with a cloth like feel.

### POLYETHELENE

A film coating applied to the non-woven polypropylene fabric offering maximum barrier properties and is ideal for any situation where a high level of fluid is anticipated.

## FEATURES


- Light weight
- Lint free
- Long sleeves
- Knitted cuffs
- Latex free
- Closed back for maximum protection
- Impervious

AVAILABLE COLOURS:



Our gowns are:

- Infective Agents tested to EN 14126.2003

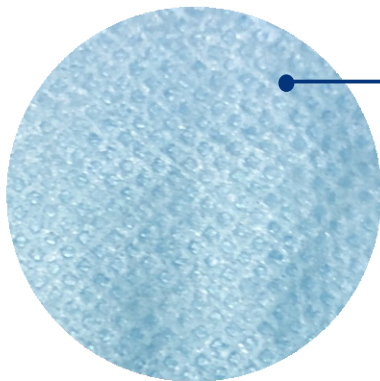

PE film coating

